




Pole Ideal Tajhiz Co.  
www.medpit.com

## Hematology Roller Mixer PITR8

Mixing and keeping the laboratory suspensions uniform

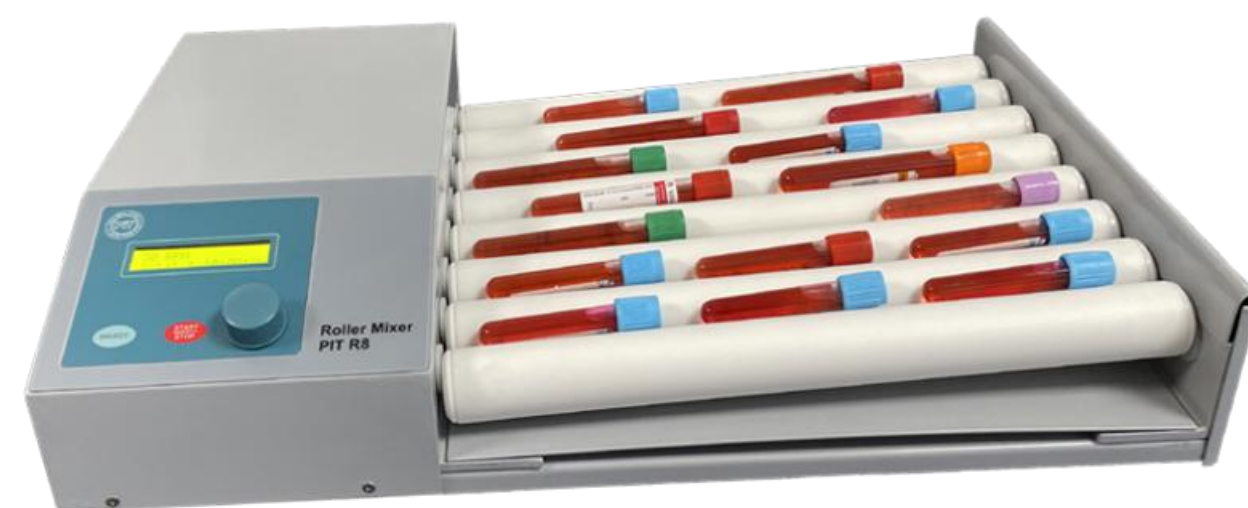

P.I.P. Hematology Mixer PITR8 has 8 rollers for mixing blood specimens, high viscosity samples, suspensions, etc. One of the functional properties of this device is its simultaneous rotational and vertical motion of the rollers. The roller dimensions and the distances between them make the instrument

usable with different tubes and cylindrical bottles. A washable protective tray is designed under the rollers in order for easy cleaning in case of sample spillage. P.I.P. Roller Mixer is resistant to most of the common laboratory acids and chemicals.

### Technical Data

Vertical motion range	~20 mm
Roller dimensions (diameter × length)	φ32 mm × 320 mm
Speed range	20-80 rpm
Input voltage	210-230 V
Frequency	50 Hz
Max. power	15 W
Display	digital
Timer domain	00:00':00" - 99:59':59"
Permissible ambient temperature	5-40 (°C)
Max. relative humidity	80 (%)
Weight (approx.)	9 kg

GS1 Code	Cat No.	Dimensions
6265738001321	PITR8	500mm × 320mm × 105mm

1<sup>st</sup> Floor, NO.12, Naghdi St., Jahantab St., Motahari Ave., Tehran, Iran

Postal Code: 1576635714

PO.Box: 15875 - 9483

✉ Info@medpip.com

🏢 Pole Ideal Pars Co.

☎ +98 9123340197

✈ @poleideal

📷 poleideal

Tel: +98 21 88545922-9

Fax: +98 21 88765561

+98 21 88767159